Medical Measuring Systems and Scales since 1840

seca gmbh & co. kg
Hammer Steindamm 3–25
22089 Hamburg · Germany
T +49 40 20 00 00 0
F +49 40 20 00 00 50
E info@seca.com

seca operates worldwide with headquarters in Germany and branches in:

- seca france
- seca united kingdom
- seca north america
- seca schweiz
- seca zhong guo
- seca nihon
- seca mexico
- seca austria
- seca polska
- seca middle east
- seca suomi
- seca américa latina
- seca asia pacific
- seca danmark
- seca benelux
- seca lietuva


and with exclusive partners in more than 110 countries.
All contact data at seca.com

17-10-05-414-100_2026-09S

Instructions for Use
seca 441 Solar Panel

1. PRODUCTS IN THIS MANUAL / INTENDED USE • PRODUITS VISÉS PAR CETTE NOTICE / DOMAINE D'UTILISATION
• PRODUCTOS EN ESTAS INSTRUCCIONES / USO PREVISTO • الغرض من الاستخدام • المنتجات في هذا الدليل

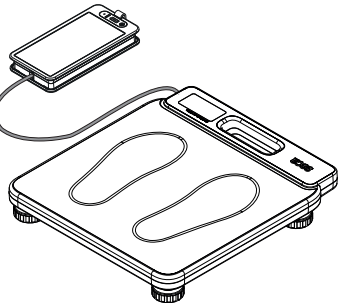The solar panel serves as a mobile power source for compatible seca scales.
For use with seca scales **seca 874 u+** and **seca 874 +** only.
Scales to be ordered separately.

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⚠ WARNING!
Personal injury, damage to device
Handling the device incorrectly or using a damaged device may lead to severe irreversible or fatal injuries.
► Do not expose device to heat sources or fire.
► Do not store device in a vehicle or in other locations subject to high temperatures.
► Do not open housing and do not repair device yourself.
► Do not continue using device if it leaks or if other damage is visible.
► Keep device out of the reach of children.
► Device is intended for use by persons in good physical and mental condition who have the necessary experience and knowledge to handle the device.

⚠ NOTICE!
Damage to device
Handling the device incorrectly can lead to damage to the device.
► Do not expose the device to shocks and do not drop it from a significant height.

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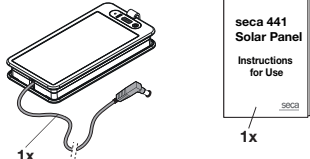
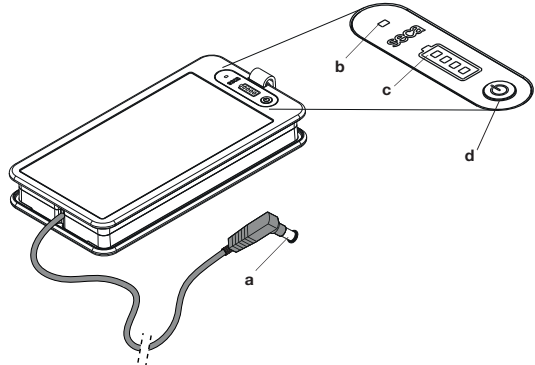
c

d

Charging LEDs	Capacity
ON OFF	0 %
ON OFF	25 %
ON OFF	50 %
ON OFF	75 %
ON OFF	100 %

Pos.	EN	FR	ES	العربية
a	Device connector	Connecteur d'alimentation	Clavija del aparato	موصل الجهاز
b	Power LED	Témoin « État de fonctionnement »	LED "estado de funcionamiento"	مؤشر الطاقة
c	Charge status	État de charge	Estado de carga	مؤشر حالة الشحن
d	On/off switch	Interrupteur marche/arrêt	Interruptor de encendido/apagado	مفتاح التشغيل/الإيقاف

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Cleaning/Nettoyage/Limpieza/التنظيف

NOTICE:

- Do not use aggressive or abrasive cleaning agents.
- Use a soft cloth dampened with mild soapsuds to wipe down the device with it.

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Disinfecting/Désinfection/Desinfección/التعقيم

NOTICE:

- Use alcohol-based disinfectant (e.g. 70 % ethanol) exclusively.
- Follow the instructions on the disinfectant.
- Moisten a soft cloth with disinfectant and wipe down the device with it.

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