
seca 952

Instructions for Use

17-10-06-660-100a_2026-09S

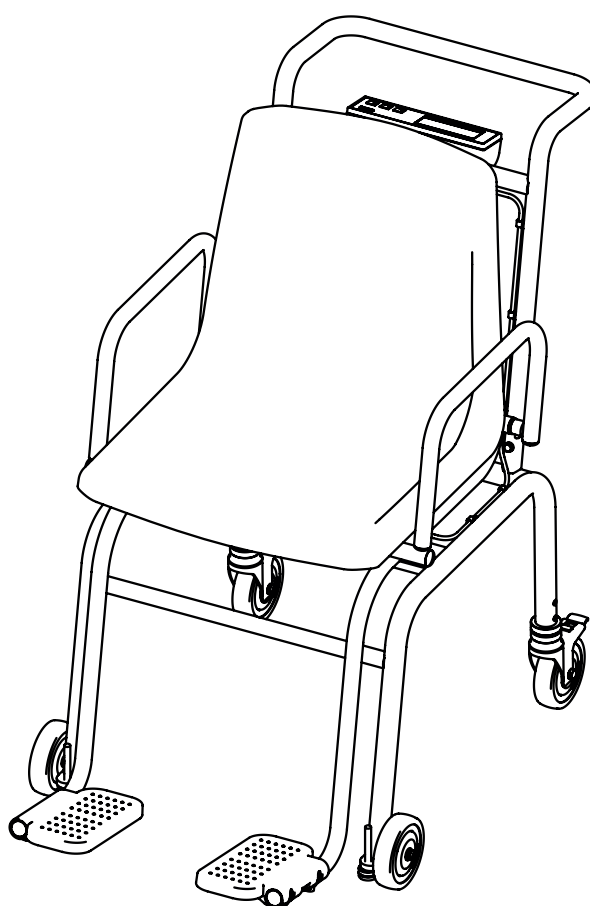Deutsch – Gebrauchsanweisung	3
English – Instructions for Use.....	22
Français – Mode d’emploi.....	42
Español – Instrucciones de uso	62
Italiano – Istruzioni per l’uso	81
Português – Manual de instruções de utilização	100
Polski – Instrukcja obsługi.....	120

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TABLE OF CONTENTS

1 About this document	23	8 Function check	37
1.1 Representation in text	23	9 Troubleshooting.....	38
1.2 Representation in diagrams	23	10 Servicing.....	39
2 Description of device	23	10.1 Servicing	39
2.1 Intended use.....	23	11 Technical data.....	39
2.2 Contraindications.....	23	11.1 General technical data.....	39
2.3 Clinical benefit	24	11.2 Weight measurement	40
2.4 Patient target group.....	24	12 Disposal	40
2.5 User qualification	24	12.1 Disposing of the device	40
2.6 Functional description.....	24	12.2 Disposing of batteries and rechargeable batteries	40
Measuring weight	24	13 Warranty	40
3 Safety precautions.....	24	14 Declaration of conformity	41
3.1 Safety information in these instructions for use	24		
3.2 Basic safety information.....	25		
Handling device	25		
Preventing injuries and infections	26		
Preventing damage to device.....	27		
Handling measuring results.....	27		
Handling packaging material.....	27		
Handling batteries and rechargeable batteries	28		
4 Overview.....	29		
4.1 Controls.....	29		
4.2 Symbols in the display	30		
4.3 Labels.....	30		
5 Starting up device.....	32		
5.1 Inserting batteries	32		
5.2 Scope of delivery	32		
6 Operation.....	33		
6.1 Setting up and aligning the device	33		
6.2 Switching the device on and off.....	34		
Switching on the device.....	34		
Switching off the device.....	34		
6.3 Weighing a patient	34		
6.4 Auto hold function	34		
6.5 Switching the weight unit.....	35		
6.6 Taring off additional weight (TARE).....	35		
7 Hygiene treatment	36		
7.1 Cleaning	36		
7.2 Disinfecting.....	36		
7.3 Sterilizing	37		

1 ABOUT THIS DOCUMENT

1.1 Representation in text

Symbol	Description
✓	Requirement for actions
►	Action
1. 2.	Actions with specified sequence
a) b)	Steps of an action with specified sequence
⇒	Result of an action
• •	First level of a list
– –	Second level of a list

1.2 Representation in diagrams

Symbol	Description
	Indicates directions of movement of the device or of device components
	Navigation path in menu displays
	• Correct action • Correct result of action
	• Incorrect action • Incorrect result of action
	End of a procedure, e.g. the installation of a part

2 DESCRIPTION OF DEVICE

2.1 Intended use

The electronic chair scale supports physicians in decision-making regarding diagnosis or therapy, based on weight or weight changes.

2.2 Contraindications

No contraindications are known.

2.3 Clinical benefit

The electronic chair scale supports physicians in decision-making regarding diagnosis or therapy, based on the measured weight.

By that, the electronic chair scale only provides an indirect clinical benefit.

The physician must perform additional examinations for an accurate diagnosis.

2.4 Patient target group

The scale is intended for persons of all ages (with exception of babies), who are able to independently sit on the chair scale without additional external support and who do not exceed the maximum weight capacity of the scale.

2.5 User qualification

The device may only be operated by persons with formal training in the healthcare sector or in medicine.

2.6 Functional description

Weight calculation is performed by a single load cell.

3 SAFETY PRECAUTIONS

3.1 Safety information in these instructions for use

DANGER!

Used to identify an extremely hazardous situation. If you fail to take note of this information, serious irreversible or fatal injuries will occur.

WARNING!

Used to identify an extremely hazardous situation. If you fail to take note of this information, serious irreversible or fatal injuries may result.

CAUTION!

Used to identify a hazardous situation. If you fail to take note of this information, minor to moderate injuries may result.

NOTICE!

Used to identify possible incorrect usage of the device. If you fail to take note of this information, you may damage the device, or the measuring results may be incorrect.

NOTE

Includes additional information about use of the device.

3.2 Basic safety information

Handling device

- ▶ Please take note of the information in these instructions for use.
- ▶ Keep the instructions for use in a safe place. The instructions for use are a component of the device and must be available at all times.
- ▶ In the interest of patient safety, you and your patients are obliged to report serious events that occur in connection with this product to the manufacturer and to the authority responsible in your country.

DANGER!

Risk of explosion

- ▶ Do not use the device in an environment in which one of the following gases has accumulated:
 - oxygen
 - flammable anesthetics
 - other flammable substances/air mixtures

CAUTION!

Patient hazard, damage to device

- ▶ Additional devices connected to medical electrical devices must provide evidence of compliance with the relevant IEC or ISO standards (e.g. IEC 60950 for data-processing devices). Furthermore, all configurations must comply with the requirements of standards for medical systems (see IEC 60601-1-1 or Section 16 of edition 3.1 of IEC 60601-1 respectively). Anyone connecting additional devices to medical electrical devices is considered a system configurator and therefore responsible for ensuring that the system complies with the requirements of standards for systems. This also applies to additional devices recommended by seca. Your attention is drawn to the fact that local laws take precedence over the above-mentioned requirements of standards. In the event of any queries, please contact your local specialist dealer or Technical Service.

CAUTION!

Patient hazard, damage to device

- ▶ Have servicing carried out regularly as described in the relevant section of this document.
- ▶ Technical modifications may not be made to the device. The device does not contain any user-serviceable parts. Only have servicing and repairs performed by an authorized seca service partner. You can find your local service partner at www.seca.com.
- ▶ Use only seca original accessories and spare parts, otherwise seca will not grant any warranty.

CAUTION!

Patient hazard, malfunction

- ▶ Keep medical electrical devices such as high-frequency surgical equipment at a minimum distance of approx. 1 meter to prevent faulty measurements.
- ▶ Keep high-frequency equipment such as cellphones at a minimum distance of approx. 1 meter to prevent faulty measurements.
- ▶ The actual transmission output of HF devices may require minimum distances of more than 1 meter. For details, go to www.seca.com.

Preventing injuries and infections

CAUTION!

Personal injury, damage to device

- ▶ Secure the scale against unintentional rolling away. Set the emergency brakes.
- ▶ Never place the scale on a sloped surface.
- ▶ Do not transport any people or objects with the scale.
- ▶ Move the scale only with released emergency brakes.

WARNING!

Injury from falling

People with limited mobility may fall when they sit down on the scale or stand back up.

- ▶ Route connector cables (if present) so that users and patient cannot trip over them.
- ▶ Make sure that the scale's emergency brakes are set.
- ▶ Make sure the footrests are pushed in (model-dependent) and folded up before the patient sits down on the scale or stands up again.
- ▶ Assist people with limited mobility when they are sitting down and standing up again.

CAUTION!

Injury from pinching

Hands and arms can get pinched between the seat and the lower frame or between the seat and the armrests.

- ▶ Keep your arms and hands above the seat when assisting the patient.
- ▶ Make sure that the patient's hands and arms are on the armrests or in their lap during weighing.

WARNING!

Risk of infection

- ▶ Before and after every measurement, wash your hands to reduce the risk of cross-contamination and nosocomial infections.
- ▶ Should the patient have contagious diseases, subject the device to a hygiene treatment as described in the relevant section of these instructions for use immediately afterwards.
- ▶ Ensure that the patient does not have any open wounds or infectious skin alterations which may come into contact with the device.
- ▶ Subject the device to a hygiene treatment at regular intervals.

Preventing damage to device

NOTICE!

Damage to device

- ▶ Ensure that fluids never get inside the device. These can destroy the electronics.
- ▶ If you are not going to use the device for an extended period of time, remove the batteries or battery packs. Only then is the device de-energized.
- ▶ Do not drop the device.
- ▶ Do not subject the device to shocks or vibrations.
- ▶ Perform a function check before each use as described in the corresponding section in this document. Do not operate the device if it is not working properly or is damaged.
- ▶ Do not place the device in direct sunlight and ensure that it is not placed in the direct proximity of a heat source. The excessive temperatures could damage the electronics.
- ▶ Avoid rapid temperature fluctuations. If the device is transported so that a temperature difference of over 20 °C occurs, the device must be left to stand for at least 2 hours before it is switched on, otherwise condensation may form; this may damage the electronics.
- ▶ Use the device only in the intended ambient conditions.
- ▶ Store the device only in the intended storage conditions.
- ▶ Use only cleaning agents and disinfectants which match the details in the section entitled "Hygiene treatment".
- ▶ For scales: Ensure that maximum capacity is not exceeded.

Handling measuring results

CAUTION!

Patient hazard

To prevent misinterpretations, measuring results for medical purposes must only be displayed and used in SI units (kilograms/grams, meters/centimeters). Some devices have the option of displaying measuring results in different units. This is purely an additional function.

- ▶ Only use the measuring results in SI units.
- ▶ You as the user take sole responsibility for the use of measuring results in non-SI units.

NOTICE!

Inconsistent measuring results

- ▶ Before you save and continue using measured values determined using this device (e.g. in seca software or in an EMR system), ensure that the measured values are plausible.
- ▶ If measured values have been sent to seca software or to an EMR system, ensure before continuing to use them that the measured values are plausible and assigned to the correct patient.

Handling packaging material

WARNING!

Risk of suffocation

Packaging materials made of plastic film (e.g. bags) present a risk of suffocation.

- ▶ Store packaging material out of the reach of children.
- ▶ If the original packaging material is no longer available, only use plastic bags with safety holes to reduce the risk of suffocation. Use recyclable materials if possible.

NOTE

Store the original packaging material for future use (e.g. returning for servicing).

Handling batteries and rechargeable batteries

WARNING!

Injury

Batteries and rechargeable batteries contain harmful substances which may explode and so escape if not handled properly.

- ▶ Do not try to recharge batteries.
- ▶ Do not expose batteries/rechargeable batteries to heat.
- ▶ Do not burn batteries/rechargeable batteries.
- ▶ If harmful substances are leaking out, avoid contact with the skin, eyes, and mucous membranes. Rinse affected areas of the body with plenty of clean water and seek medical help at once.

NOTICE!

Damage to device and malfunction due to incorrect handling

- ▶ Use only the type of battery/rechargeable battery quoted in this document.
- ▶ Always replace all batteries/rechargeable batteries at the same time.
- ▶ Do not short-circuit batteries/rechargeable batteries.
- ▶ If you do not use the device for an extended period of time, remove the batteries/rechargeable batteries. This prevents acid from leaking into the device.
- ▶ If harmful substances have penetrated the device, do not continue using it. Have the device checked by an authorized seca service partner and repaired if necessary.


4 OVERVIEW

4.1 Controls

Item	Control	Function
1	Armrests	Can be folded back
2	Battery compartment	Takes 6 AA-type mignon batteries (AA/1.5 V)
3	Footrests	Can be folded up
4	Brake levers	Used to secure the rear rollers
5		Used for switching the device on and off
6		kg/lbs/sts key - Press briefly: Switches between kg and lbs - Keep depressed: Switches between lbs and sts
7		hold/tare key - Press briefly: Activate Hold function - Press and hold: Activate Tare function

4.2 Symbols in the display

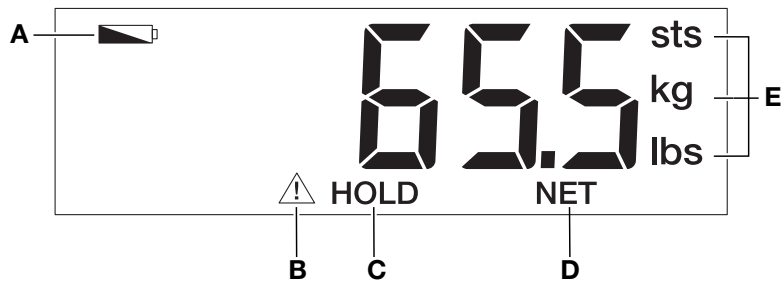Item	Symbol	Meaning
A		Batteries are low
B		Non-verifiable function active (for technical reasons, symbol is also displayed on non-verified scales)
C	HOLD	Hold function active
D	NET	Tare function active
E	kg/lbs/sts	Unit in which weight is displayed

4.3 Labels

Markings on the device and on the type plate	
Text/symbol	Meaning
	Name and address of manufacturer, date of manufacture
UDI	Unique Device Identification (product identification number in accordance with Regulation (EU) 2017/745)
REF	Article number
SN	Serial number
ProdID	seca product identification number
Mat.No.	Variant number
	Follow instructions for use
	Medical electrical device, type B
IP20	Type of protection in accordance with IEC 60529: • Protection against ingress of solid foreign objects with a diameter of over 12.5 mm • Protection against access with fingers • Protection against dripping water
d	Value in units of mass States the difference between two consecutive display values
MD	Medical device in accordance with Regulation (EU) 2017/745

Markings on the device and on the type plate	
Text/symbol	Meaning
	Device complies with EU directives 0123 : Notified Body for Medical Devices
	Device meets the requirements of the USA and Canada. Certified and tested by a licensing laboratory (NRTL) of TÜV SÜD Product Services GmbH.
	Device complies with United Kingdom directives XXXX : Notified Body for Medical Devices of the United Kingdom
	Importer/representative in the United Kingdom: seca Ltd 40 Barn Street B5 5QB Birmingham United Kingdom
	Importer/representative in Switzerland: seca ag (schweiz) Medizinische Waagen und Messsysteme Schönmatt Str. 2 CH-4153 REINACH
	Do not dispose of device in household waste

Markings on the packaging	
Symbol	Meaning
	Protect from moisture
	Arrows indicate top of product Transport and store in an upright position
	Fragile, do not throw or drop
	Permitted min. and max. temperature for transport and storage
	Permitted min. and max. humidity for transport and storage
	Permitted min. and max. air pressure for transport and storage
	Open packaging here
	Packaging material can be disposed of through recycling programs

5 STARTING UP DEVICE

5.1 Inserting batteries

NOTICE!

Malfunction as a result of errors in installation

If cables are strained during installation, e.g. with sharp bends or kinked connectors, this may result in faulty displays and failure of the display.

- Route all cables to prevent sharp bends and kinked connectors.
- Provide strain relief by routing all cables in the relevant holders.

✓ Six batteries (type AA/1.5 V)

1. Open the battery compartment.
2. Connect the cable to the battery holder.


NOTICE!

Damage to device


The device may overheat and become damaged.

- Ensure the correct polarity of the batteries (markings in the battery compartment).

3. Place the batteries in the battery compartment.
4. Close the battery compartment.

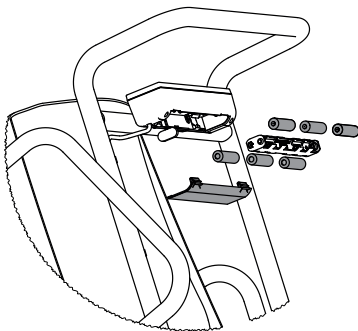5.2 Scope of delivery

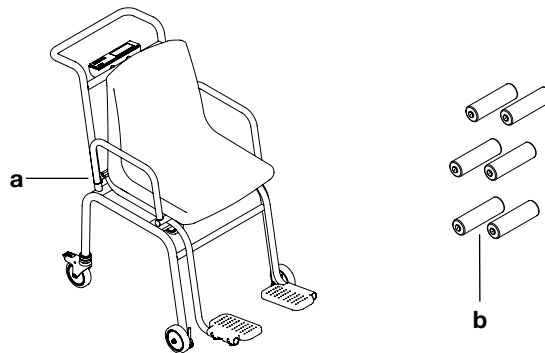Item	Component	Pcs.
a	Chair scale (fully assembled)	1
b	Batteries (type AA/1.5 V)	6
Not shown.	Instructions for use	1

6 OPERATION

WARNING!

Personal injury, faulty measurement

- Before using the device each time, perform a function check as described in the corresponding section of these instructions for use.

WARNING!

Injury from falling

People with limited mobility may fall when they sit down on the scale or stand back up.

- Route connector cables (if present) so that users and patient cannot trip over them.
- Make sure that the scale's emergency brakes are set.
- Make sure the footrests are pushed in (model-dependent) and folded up before the patient sits down on the scale or stands up again.
- Assist people with limited mobility when they are sitting down and standing up again.

CAUTION!

Injury from pinching

Hands and arms can get pinched between the seat and the lower frame or between the seat and the armrests.

- Keep your arms and hands above the seat when assisting the patient.
- Make sure that the patient's hands and arms are on the armrests or in their lap during weighing.

CAUTION!


Faulty measurement as a result of force shunt

If accessories or articles of clothing get caught between the seat and the frame, the weight will not be measured correctly.

If the patient's feet are not on the footrests, the weight will not be measured correctly.

- Make sure that there are no objects between the seat and frame.
- Make sure that the armrests are folded down during weighing.
- Make sure that the patient's feet are on the footrests during weighing.


6.1 Setting up and aligning the device

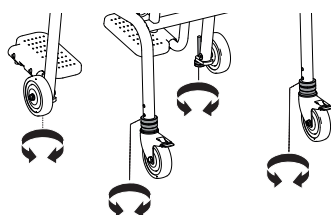1. Place the scale on a firm, level surface.
2. Loosen the knurled head screws on the foot screws and on the rear rollers (rear rollers: loosen the upper knurled head screws).
3. Align the scale:
 - Rear rollers: turn the lower knurled head screws
 - Front rollers: turn the foot screws
4. Make sure all four of the scale's rollers are in contact with the ground.
5. Secure the adjustment on all four rollers by tightening the knurled head screws (rear rollers: tighten the upper knurled head screws)

6.2 Switching the device on and off

Switching on the device

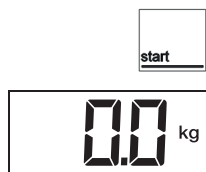1. Ensure that the device has no load.
2. Press the **START** key.
 - ⇒ All the elements of the display are shown briefly.
 - ⇒ **SECA** is displayed.
 - ⇒ The device is ready for operation when **0.0** is displayed.

Switching off the device


- Press the **START** key.

NOTE

In battery mode, the device remains on standby for 3 minutes. If there is no input during the standby time, the device switches off automatically.

6.3 Weighing a patient

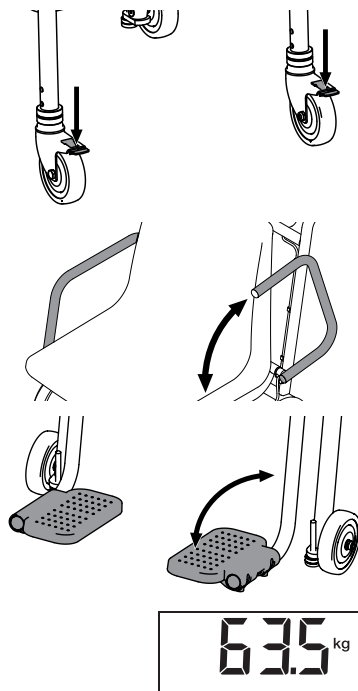- ✓ The scale has no load.
 - ✓ The scale is switched on.
1. Set the emergency brakes of the rear casters.
 2. Fold the armrest toward the back in order to make it easier for the patient to sit down.
 3. Ask the patient to sit down on the scale.
 4. Fold the footrests down.
 5. Ask the patient to put their feet on the footrests.
 6. Fold the armrests down.
 7. Read off the measuring result.

6.4 Auto hold function

The **Hold** function is activated automatically at weights of over 1.5 kg. The display flashes until a stable weight has been measured. The display then freezes until the next weighing operation. The message **HOLD** is displayed.

- To take another measurement without first unloading the scale, press the **hold/tare** key.

6.5 Switching the weight unit

You can select the unit in which weight is to be displayed; kilograms (**kg**), pounds (**lbs**) or stones (**sts**).

CAUTION!


Patient hazard

To prevent misinterpretations, measuring results for medical purposes must only be displayed and used in SI units (kilograms/grams, meters/centimeters).

- ▶ Only use the measuring results in SI units.
- ▶ You as the user take sole responsibility for the use of measuring results in non-SI units.


✓ The scale has no load.

1. Press the **START** key.

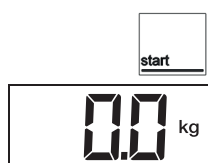⇒ The current weight unit will be displayed (in this case **kg**).

2. To switch the weight value between kilograms, pounds, and stones, press the **kg/lbs/sts** key repeatedly until the desired unit appears on the display.

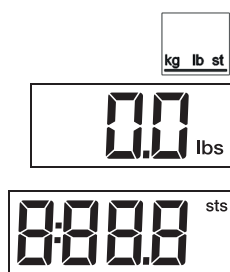6.6 Taring off additional weight (TARE)

Using the **TARE** function, you can prevent additional weight (e.g. a hand towel or a support on the seat) from influencing the weighing result.

CAUTION!

Faulty measurement as a result of force shunt

If accessories or articles of clothing get caught between the seat and the frame, the weight will not be measured correctly.

If the patient's feet are not on the footrests, the weight will not be measured correctly.


- ▶ Make sure that there are no objects between the seat and frame.
- ▶ Make sure that the armrests are folded down during weighing.
- ▶ Make sure that the patient's feet are on the footrests during weighing.

NOTE

The maximum weight which can be displayed is reduced by the weight of the objects already placed on the scale.

✓ The scale has no load.

1. Press the **START** key.
2. Place the additional weight on the scale.
3. Press and hold the **hold/tare** key.
 - ⇒ The display is set to zero.
 - ⇒ The message **NET** is displayed.

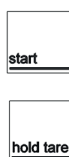**NOTE**

The patient's weight can only be measured when they are no longer in contact with the ground.

- Make sure that the patient's feet are on the footrests during weighing.

4. Ask the patient to sit down on the scale.
5. Ask the patient to put their feet on the footrests.
 - ⇒ The patient's weight is displayed.
 - ⇒ The additional weight is deducted automatically.
6. Read off the measuring result.

**NOTE**

Now, you can perform as many measurements as you like with this additional weight.

7. To deactivate the **Tare** function, briefly press the **hold/tare** key.
 - ⇒ The **NET** message is no longer displayed.

NOTE

When you switch the device off, the **Tare** function is deactivated automatically.

7 HYGIENE TREATMENT

**WARNING!****Electric shock**

The device is not de-energized when the on/off button is pressed and the display goes out. Use of fluids on the device may cause an electric shock.

- Before each hygiene treatment, ensure that the device is switched off.
- Take the batteries/battery packs out of the device before every hygiene treatment.
- Ensure that no fluids penetrate the device.

**WARNING!****Risk of infection**

- Subject the device to a hygiene treatment at regular intervals as described in this section.

NOTICE!**Damage to device**

Unsuitable cleaning agents and disinfectants may damage the sensitive surfaces of the device and impair its operability.

- Do not use aggressive or abrasive cleaning agents.
- Do not use organic solvents (e.g. white spirit or petroleum spirit).

7.1 Cleaning

- If required, moisten a soft cloth with a mild soap solution and wipe the device over with it.

7.2 Disinfecting

1. Disinfect the device at regular intervals with a disinfectant suitable for sensitive surfaces and acrylic glass (e.g. 70 % ethanol).
2. Follow the instructions for use of the disinfectant.

3. Disinfect the device:

- ▶ Moisten a soft cloth with disinfectant and wipe down the device with it.
- ▶ Comply with the intervals, see table:

Interval	Component
Before and after every measurement with direct skin contact	• Seat • Armrests • Footrests
As required	• Display • Keyboard • Handle bar

7.3 Sterilizing

The device must not be sterilized.

8 FUNCTION CHECK

CAUTION!

Personal injury

If you find faults or deviations during the function check which you are unable to remedy with the aid of the "Troubleshooting" section in this document, you must not use the device.

- ▶ Have the device repaired by seca Service or an authorized service partner.
- ▶ Follow the section entitled "Servicing" in this document.

- ▶ Perform a function check before each use.

A complete function check includes:

- A visual inspection for mechanical damage
- A test of the alignment of the device
- A visual and function check of the display elements
- A function check of all the controls shown in the section entitled "Overview"
- A function check of the optional accessories

If you find faults or deviations during the function check, initially try to remedy the fault with the aid of the "Troubleshooting" section in this document.

9 TROUBLESHOOTING

Fault	Cause	Remedy
No weight is displayed with a load on the scale.	The device has no power supply.	- Switch on device. - Insert batteries.
 is displayed.	Battery voltage is dropping.	Insert new batteries as soon as possible
bAtt is displayed.	The batteries have discharged.	Insert new batteries
0.0 is not displayed before weighing takes place.	A load was already on the device before it was switched on.	- Remove the load from the device. - Switch off the device, then switch it back on again.
STOP is displayed.	Maximum capacity has been exceeded.	Remove the load from the device.
The device no longer reacts when keys are pressed.	The keypad is defective.	Inform seca Service.
One display segment is displayed continuously or not at all.	The display is defective.	Inform seca Service.
Implausible measuring results.	Patient moved too much during the measurement.	Take the load off the device, wait until 0.0 is displayed and repeat the measurement procedure.
The weight display flashes.	Consistent weight value not yet detected.	- One of the following functions is active: Hold, Tare. Wait until the device has detected a consistent weight value. The display stops flashing. - Take the load off the device, wait until 0.0 is displayed and measure again.
----- is displayed before weighing takes place.	A load was already on the device before it was switched on.	- Remove the load from the device. - Switch off the device, then switch it back on again.
Er:[No.]:11 is displayed.	The device has too high a load or too high a load in one corner.	- Remove the load from the device or distribute the weight more evenly. - Switch off the device, then switch it back on again. - If the problem persists, inform seca Service.
Er:[No.]:12 is displayed.	The device has been switched on with too high a load.	- Remove the load from the device. - Switch off the device, then switch it back on again.
Er:[No.]:16 is displayed.	The device was caused to oscillate, the zero point could not be determined.	Switch off the device, then switch it back on again.
Er:[No.]:40 is displayed.	Error in the configuration memory	Inform seca Service
Er:[No.]:50 is displayed.	Error in program memory	Inform seca Service

10 SERVICING

10.1 Servicing

The product needs to be set up carefully and serviced regularly. Depending on how frequently the product is used, seca recommends servicing at intervals of 3 to 5 years.

CAUTION!

Faulty measurements as a result of poor servicing

- ▶ Have servicing and repairs carried out exclusively by seca Service or an authorized service partner.
- ▶ You can find your local service partner at www.seca.com.

11 TECHNICAL DATA

11.1 General technical data

Parameter	Value
Dimensions • Depth (footrests pulled out) • Depth (footrests pushed in) • Width • Height • Seat height	978 mm 808 mm 563 mm 897 mm 528 mm
Weight	approx. 25 kg
Height of digits, display	25 mm
Ambient conditions, operation • Temperature • Air pressure • Humidity	+10 °C to +40 °C / +50 °F to 104 °F 700 – 1060 hPa 30 % – 80 % no condensation
Ambient conditions, storage/transport • Temperature • Air pressure • Humidity	-10 °C to +65 °C / +14 °F to 149 °F 700 – 1060 hPa 0 % – 95 % no condensation
Medical device in accordance with regulation (EU) 2017/745	Class I with measuring function
IEC 60601-1: Medical electrical device, Type B	
Type of protection in accordance with IEC 60529	IP20
Duty cycle	Continuous duty
Service life, from commissioning ^a	at least 10 years
Application parts according to IEC 60601-1	Seat, armrests, footrests
Batteries • Supply voltage • Battery type • Maximum running time	9 V 6 x type AA 1.5 V approx. 10000 weighing operations

^a Minimum service life when used as intended and serviced in a timely manner. seca products usually have a significantly longer service life.

11.2 Weight measurement

Non-verified model	
Maximum capacity	200 kg/440 lb/32 st
Minimum capacity	1 kg
Graduation	100 g/0.2 lb
Tare range	200 kg (subtractive)
Accuracy	
• < 100 kg	± 100 g
• 100 kg to 200 kg	± 0.1 %
• < 200 lb	± 0.2 lb
• 200 lb to 440 lb	± 0.1 %
• < 14 st	± 0.2 lb
• 14 st to 32 st	± 0.1 %
At ambient temperature of 20° C and humidity of 45 %	

12 DISPOSAL

12.1 Disposing of the device

Do not dispose of the device in your household waste. The device must be properly disposed of as electronic scrap. Follow your respective national regulations. For more information, please contact seca Service or an authorized service partner.

12.2 Disposing of batteries and rechargeable batteries

Do not dispose of used batteries and rechargeable batteries in household waste, regardless of whether they contain harmful substances or not. As a consumer, you are legally obliged to dispose of batteries and rechargeable batteries via local authority collection points or trade collection points. Only dispose of batteries and rechargeable batteries once they are completely discharged.

13 WARRANTY

There is a two-year warranty period from delivery for defects attributable to poor materials or workmanship. All movable parts, e.g. batteries, cables, power supply units, rechargeable batteries etc. are exempt. Defects which come under the warranty will be repaired for the customer free of charge against proof of purchase. Additional claims cannot be considered. Costs of transport to and from seca are the responsibility of the customer if the device is located somewhere other than the customer's headquarters. In the event of transport damage, claims under warranty can only be made if the complete original packaging was used for transport and the device was secured and fastened in it according to its originally packaged condition. You should therefore keep all packaging parts.

The warranty will be voided if the device is opened by persons not expressly authorized by seca to do so.

In the event of a warranty issue, please contact your local seca office or the dealer from whom you ordered the product.

14 DECLARATION OF CONFORMITY

seca gmbh & co. kg hereby declares that the product complies with the terms of the applicable European directives and regulations. The unabridged declaration of conformity can be found at www.seca.com.

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