seca 777 series



Time to say good buy!



seca 777 series The price-performance


champion of measuring and weighing

You don't have to compromise on precision or quality when looking for a cost-efficient scale. The seca 777 provides the perfect solution with its oversized low-profile, tiltproof solid metal platform (no plastic), builtin measuring rod with 4 – 90" measuring range and optional integrated handrail.

The design sets standards for precise, quick and secure weighing and measuring. With its high capacity of 550 lbs, the scale is optimal for accommodating all patient types. Its large, rubber coated castors also make it convenient for staff to transport the scale when needed.

If you are looking to upgrade to a digital measuring rod or EMR integration, the seca 787 and 797 are the perfect options for you.

_ow-pro	file solid metal platform with tilt-proof design
Measurir	ng rod (4 – 90 inch / 10 – 230 cm)
_arge ru	ober coated castors for excellent movability
Display s	swappable to the left or right side
Optional	integrated handrail accessory
Optional	panda facilitates measuring of children
Damping	system and Auto-HOLD function
Digital m	easuring rod (from 48 inch / 122 cm up)
Two-line	display shows height, weight and BMI
JSB-por	t and holder for scanner*
EMR val	dated with upgradeable LAN/Wi-Fi module
EMR val	idated with built-in Wi-Fi module

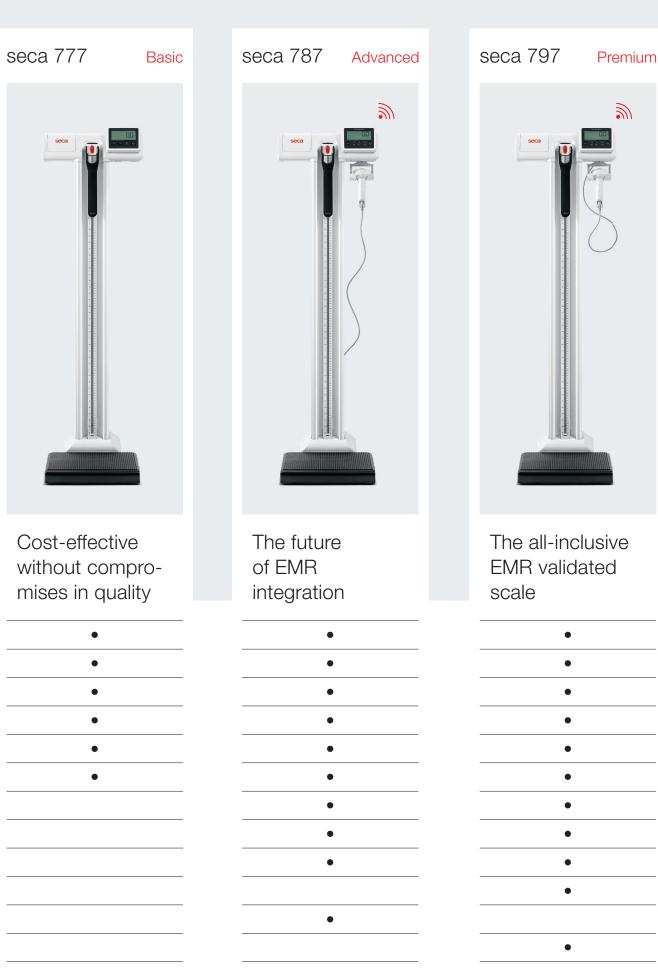





* Scanner from a third-party vendor





seca integration solutions A clear improvement to your workflow

The seca 787 and seca 797 can be integrated into your information system via the existing network. At the push of a button, measured values can be sent directly to the electronic medical record (EMR), where they are immediately available to all users.

Connected with an optional USB barcode scanner, both patient and user are identified and the measurement results are always correctly assigned. This saves time, increases process reliability and eliminates transcription errors.

Get further information on seca integration solutions at: integration.seca.com



Technical data	seca 777	seca 787	seca 797	
Functions	TARE, HOLD, lbs/kg switch-over lock, BMI	TARE, HOLD, Auto-HOLD, Damping (motion filter), lbs/kg switch-over lock, automatic BMI, auto-off	TARE, HOLD, Auto-HOLD, Damping (motion filter), kg/lbs switch-over lock, automatic BMI, USB-port and holder for scanner	
Data transmission	USB interface for customer specific solution	LAN or Wi-Fi via interface module seca 452 (optional)	Integrated Wi-Fi module	
Power supply	Batteries, power adapter (optional)	Batteries, power adapter (optional)	Power adapter, batteries (optional)	
Accessories	Handrail seca 477, Panda seca 459, Power adapter 401	Handrail seca 477, Panda seca 459, Power adapter 401	Handrail seca 477, Panda seca 459	
Weight	28.0 lbs / 12.7 kg	28.4 lbs / 12.9 kg	28.9 lbs / 13.1 kg	
Measuring range (total)	4 – 90 inch / 10 – 230 cm	4 – 90 inch / 10 – 230 cm	4 – 90 inch / 10 – 230 cm	
Measuring range (digital)		48 – 90 inch / 122 – 230 cm	48 – 90 inch / 122 – 230 cm	
Capacity	550 lbs / 250 kg			
Increments	0.2 lbs / 100 g			
Dimensions (WxHxD)	16.1 x 53.4 x 22.7 inch / 41.0 x 135.6 x 57.6 cm			
Platform (WxHxD)	13.8 x 2.3 x 16.5 inch / 35.0 x 5.9 x 41.8 cm			
Electromagnetic radiation resistance	60601-1-2 compliant			

201377704020 seca north america 13601 Benson Ave Chino, CA 91710

USA

+1 800 542 7322 Т F +1 888 705 7397 E info.us@seca.com

Precision for health seca.com/777