# SCALING UP

# WEIGHING& MEASURING REDEFINED

seca SCALE-UP LINE ——





Play video to see the scale in action!

seca Mobile Handrail Scale

# Built for Durability Designed for Safety Ready for High Traffic

The seca Scale-up Line is purpose-built to endure the rigors of daily medical work with resilience. Its heavy-duty, robust design ensures consistent accuracy even after years of relentless service, providing healthcare teams with trusted reliability.



#### MOBILITY MEETS DURABILITY

The optimized mobile and easy-to-move model of the seca Scaleup Line incorporates a digital height rod and optional rechargeable battery pack. Digital height measuring results appear in the display automatically with millimeter precision.



rechargeable battery pack provides more continuous weighing



Obtain weight and digital height conveniently at each patient room, bedside. or wherever needed.

The rock-solid steel handrail enhances stability to reduce fall



### DESIGNED TO THRIVE

- Ultrasound height measuring (1)
- (2)Sturdy aluminum column
- (3) Customizable touch-display
- (3) Voice-guided automation
- (4)Robust steel handrail
- (5) Patient and user identification\*
- (6) Optional rechargable battery pack
- (7)Heavy-duty wheels
- (8) Break-resistant Bearclaw glass
- (8) Easy to clean surfaces
- Capacity of 800 lb / 360 kg
- Fine 0.1 oz / 50 g graduation

Optional barcode scanner mount seca 463, and barcode scanner from a thirdparty vendor

## in Demanding Healthcare Environments

The seca Scale-up Line redefines medical scales' standards, delivering precision and versatility that will reshape your expectations. Its innovative features and unmatched durability set a new benchmark in the industry, making it the go-to choice for healthcare professionals.

The optional panda motif makes taking children's measurements lots of fun

#### One Scale to Master all Weighing Challenges

The seca Scale-up Line is your versatile solution, perfect for a wide range of patients and healthcare scenarios. Whether you're working with pediatric or geriatric patients, need mobile functionality, or require seamless EMR integration, this scale was designed to excel in every environment.

animal sticker accessory.



Quick digital eight and weight measurement during triage.







Bedside support for convenient weighin

## Seamless Direct **EMR-validated Integration**

Enhance efficiency, eliminate transcription errors, and foster a zero-incident environment with the seca Scale-up Line. Seamlessly integrate patient measurements into your EMR\* System. Imagine the relief for staff as patient weight and height data is automatically captured. Contact us for more information on our integration solutions.





82.50. 176.20. - 36.4



HL7

MEDITECH

\* EMR integration features are functional exclusively when the device is connected to the AC Adapter and are not applicable when powered by a mobile battery pack.



The seca ID-Display shows patient and user IDs and transfers measurement data with the push of a button.



### TAILOR YOUR SOLUTION

Find the per	fect Scale-up	Line model for	your needs
--------------	---------------	----------------	------------

eeds				
	Features	Handrail scale with ID-Display <sup>™</sup>	Mobile handrail scale with digital measuring rod and ID-Display <sup>™</sup>	Ultrasonic me station with ID
	+ Flat scale seca 654 Capacity: 800 lbs / 360 kg or 360 kg Graduation: 0.1 lb / 50 g or 50 g	•	•	•
	+ Handrail seca 455	•	•	
	<ul> <li>Digital height rod seca 254:</li> <li>4 - 90 inch / 10 - 228 cm</li> <li>(digital measuring: 48 - 90 inch / 122 - 228 cm)</li> </ul>		•	
	+ Ultrasonic height measurement seca 257: without handrail: 24 – 87 inch / 60 – 220 cm with handrail: 40 – 87 inch / 100 – 220 cm			•
	+ Rechargable battery pack seca 453 Can be retrofitted to Scale-up Line devices with serial number 10000001350445 or higher	•	•	
	+ Animal sticker seca 487			٠
	+ Panda seca 459 Not in combination with battery pack	•	•	
	+ EMR integration The features are functional exclusively when the device is connected to the AC Adapter and are not applicable when powered by a mobile battery pack	•	•	•

ON SHM INT NN

7

Ø,

ON SHM ITT NN

seca headquarters (resonsible for CIS & Africa)

+49 40 20 00 00 0 info@seca.com seca gmbh & co. kg Hamburg, Germany


seca middle east +971 50 410 2319 info.ae@seca.com

Dubai United Arab Emirates

seca asia pacific +60 173 711 973 info.aspac@seca.com Kuala Lumpur Malaysia

seca.com/scaleup

Order code

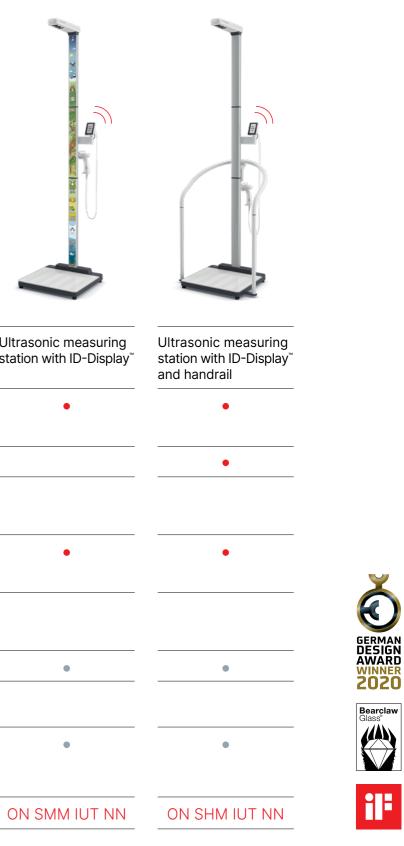

nages contain the optional barcode scanner mount seca 463 and barcode scanners from third-party suppliers.

200811