

seca mVSA 535

seca
Precision for health

Spot-check Vital Signs Monitor

We make taking vital signs EASY



seca mVSA

Where simplicity meets accuracy

We know that the working environment for medical staff is complex. This is precisely why we developed the seca mVSA, a spot-check monitor that simplifies and quickens the routine task of capturing patient vitals (blood pressure, temperature, pulse rate, SpO₂) in less than 60 seconds.

- + **EASY, quick and accurate**
- + **Configurable seca EQ BP® average mode (BP30)**
- + **seca scale connectivity through seca 360° proximity**
- + **EMR-validated and certified**
- + **Incredible price-point**
- + **World-class service program**



Optional mobile stand seca 475.

Component options



Blood pressure cuffs
(sizes XS-XXL)



Temperature probe
(oral/axillary)



Temperature probe
(rectal)



Masimo RD SET® SpO₂
finger clip



All measurements at a glance on the easy to operate 7" touchscreen monitor.



Long lasting lithium-ion battery enables up to 150 measurements on one single charge.



Store up to 70,000 vital signs measurements in the large internal memory for later processing.



Durable BP cuffs enhanced with metal plate and cuff stitching around the exterior.

EASY, quick and accurate

The seca mVSA is easy-to-use and self-explanatory with little to no training required for nursing staff. Just turn it on, start the measurement and within 60 seconds blood pressure, temperature, pulse rate and SpO₂ are shown at a glance.




Record **blood pressure** with the unique inflation or deflation mode with a push of only one button in less than 20 seconds.



The **temperature** measurement starts automatically when the probe is used (oral/axillary or rectal).



The patient's **pulse rate** is always measured simultaneously while using blood pressure measurement or pulse oximeter.

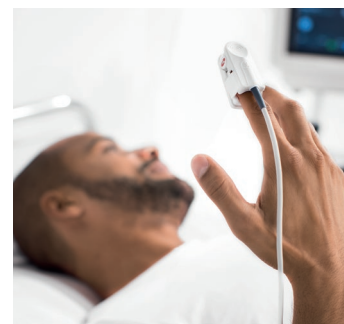

Measure **blood oxygen levels** with a Masimo RD SET® sensor.

seca EQ BP™

True screening of hypertension with configurable seca EasyQuick Blood Pressure (EQ BP™) average mode: take 6 readings automatically in 30 minutes. A single reading can be taken in under 20 seconds with the comfortable and unique blood pressure inflation mode.



Blood pressure can be determined quickly and gently using the inflation measurement method.



Avoid the “white coat effect” and capture blood pressure accurately with configurable seca EQ BP™ averaging.

Seamless scale connectivity

The seca mVSA can gather height and weight measurements from EMR-validated seca scales and measuring stations, thereby improving and standardizing base parameter acquisition.



EMR-validated and certified

Eliminate transcription errors and increase patient safety with direct transmission of measurements into the electronic medical record (EMR) via HL7. Similar integration solutions are provided for a wide range of EMR systems:

Epic **athenahealth**
eClinicalWorks

Future safe with seca service

- + Advanced replacement service can immediately provide you with an equivalently configured unit to ensure minimal downtime
- + Priority access to latest seca software updates and upgrades
- + Optimize your investment with a convenient one-, three- or five-year contract
- + Outstanding Engineering Technicians dedicated to our nationwide customers. Operating times: 9-7 EST, 8-6 CST, 6-4 PST

Toll-free service hotline: +1 800 542 7322



Configurations



Spot-check Vital Signs Monitor

- + **NiBP**
1 BP-Cuff size L
- + **TEMP (oral/axillary)**
1 Covidien Filac™ 3000
Thermometer with
oral/axillary probe
2 packages á 20 pcs
Covidien Filac™ Probe
Covers

535 1010 131



Spot-check Vital Signs Monitor

- + **NiBP**
1 BP-Cuff size L
- + **TEMP (rectal)**
1 Covidien Filac™ 3000
Thermometer with
rectal probe
2 packages á 20 pcs
Covidien Filac™ Probe
Covers

535 1020 131



Spot-check Vital Signs Monitor

- + **NiBP**
1 BP-Cuff size L
- + **SpO₂**
1 Masimo RD SET® adult
finger sensor
1 patient cable

535 1200 131



Spot-check Vital Signs Monitor

- + **NiBP**
1 BP-Cuff size L
- + **TEMP (oral/axillary)**
1 Covidien Filac™ 3000
Thermometer with
oral/axillary probe
2 packages á 20 pcs
Covidien Filac™ Probe
Covers
- + **SpO₂**
1 Masimo RD SET® adult
finger sensor
1 patient cable

535 1210 131



Spot-check Vital Signs Monitor

- + **NiBP**
1 BP-Cuff size L
- + **TEMP (rectal)**
1 Covidien Filac™ 3000
Thermometer with
rectal probe
2 packages á 20 pcs
Covidien Filac™ Probe
Covers
- + **SpO₂**
1 Masimo RD SET® adult
finger sensor
1 patient cable

535 1220 131

Accessories



Blood pressure cuffs

Size XS 5.1 – 7.9" / 13 – 20 cm	Size S 6.7 – 10.2" / 17 – 26 cm	Size M 9.4 – 12.6" / 24 – 32 cm
---------------------------------------	---------------------------------------	---------------------------------------

490 0024 001 490 0001 001 490 0002 001

Size L 12.6 – 16.5" / 32 – 42 cm	Size XL 15.0 – 18.1" / 38 – 46 cm	Size XXL 16.5 – 21.3" / 42 – 54 cm
--	---	--

490 0003 001 490 0004 001 490 0029 001



Extension tubes

Length:
51.2" / 130 cm

490 0033 001



Filac™ probe covers

Oral/axillary, rectal
2,000 pieces
(100 x 20)

490 0150 001



Mobile stand

Height adjustable:
27.4 – 46.6" /
69.5 – 118.5 cm

475 0535 009



Carrying bag

Dimensions (W x H x D):
16.1 x 12.2 x 11.6" /
41 x 31 x 29.5 cm

432 0000 009