
# seca 334 Mobile digital baby scale



- With Breast-Milk-Intake-Function (BMIF).
- Handy and light for mobile medical use.
- Comfortable cradle.
- Large, easy to read LCD display.
- Reliable results due to automatic damping system.
- Optional measuring rod seca 232.

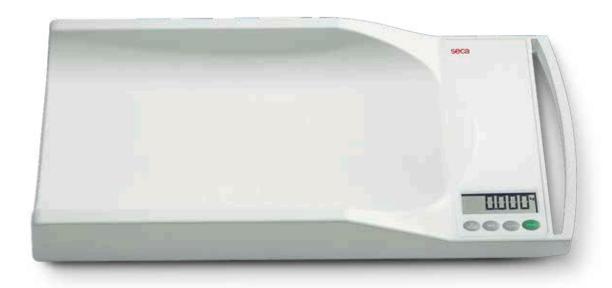

### seca 334:

# The mobile baby scale with special damping system.

odern technology can only be of use if it passes the tests of everyday life. The digital baby scale seca 334 has all that is necessary to make daily routine in the surgery and hospital so much easier as this scale is not only easy to use but also mobile and durable, it meets all the prerequisites of an ideal scale for pediatrics.

Handy and lightweight for mobile medical use.

Weighing only 6.2 lbs, this scale is a real lightweight in its class. The integrated side handle makes it easy to transport. Another practical feature is the flat and thus space-saving design, making it easy to store the scale away in the drawer of a baby's changing table. It can be packed away safely for transport in the optional transport case seca 428.



Easy to operate with generously sized cradle and easy to read LCD display.

The baby feels at ease and safe in the extra large, comfortable weighing tray of the seca 334. Larger babies or small children can be weighed in a sitting position. The weight is shown on the brightly illuminated LCD display and the automatic switch-off ensures powersaving operation. With the Breast-Milk-Intake-Function (BMIF), one can determine how much milk the baby has drank when breast fed. At the touch of a button, the weight of the baby determined before feeding is deducted from the weight determined after feeding.



Reliable results due to the automatic damping system.

Absolute precision in all situations: fitted with highly sensitive weighing technology and the special seca damping system, exact results are provided with the seca 334. Even when babies are very restless.

Optional measuring rod seca 232.

When fitted with the optional measuring rod seca 232, the baby scale seca 334 becomes a time-saving, combined measuring and weighing station.







The optional measuring rod permits weighing and measuring in one step.

.

# seca **334**

## Technical data

- Capacity: 44 lbs / 20 kg
- Graduation: 0.5 oz / 5 g < 10 kg > 10 g
- Dimensions (WxHxD): 25.1 x 4.1 x 11.8" / 638 x 105 x 300 mm
- Dimensions/weighing platform (WxHxD):
   23.4 x 2 x 10" / 595 x 50 x 255 mm
- Weight: 6.2 lbs / 2.8 kg
- Power supply:
- batteries / power adapter optional
- Functions: TARE, BMIF, auto-HOLD, lbs/kg switch-over, automatic switch-off
- Optional: measuring rod seca 232, head and foot positioner seca 419, switch mode power adapter seca 400, carrying case seca 428