seca **878**

Flat scales for mobile use



- Mother/child-function for trouble-free weighing of small children.
- Large leveling bases for increased stability.
- Low device weight and battery power for mobile use.
- Double display for readout from different positions.
- Two push buttons for effortless operation with the foot.



₩

seca 878:

The scale that thinks of people.

The seca 878 can do everything that a scale in mobile use should do and then some. It features extremely robust construction and easy-to-use design with long-life batteries and precision unaffected by humidity levels. The scale is so simple that you need only the tip of your toes to operate it.



Great stability.

The seca 878 is very light at only 4.1 kg, but also very robust. Equipped with a steel frame, the scale also boasts a top-quality ribbed surface, extra-strength two-component ABS plastic and Plexiglas[®]. Small children and adults up to 200 kg can be weighed accurately. Especially large leveling bases keep the scale stable on any type of flooring.


Simple operation.

When conducting one examination after the other, there's no need to bend over repeatedly, thanks to the two push buttons on the seca 878 that can be turned on and off with a toe tap. Yet another labor-saving feature is the mother/child-function that weighs a child held in its mother's arms. Transport of the scale is made easy by the integrated handle and very low device weight of 4.1 kg.

Intelligent details.

Special mention should be made of the double display that lets the patient and medical personnel read the results from two different directions at the same time. Thanks to a protected toggle switch and a manual power switch, which prevents unintended activation during transport, the function of the scale is always ensured even under difficult conditions.

Both doctor and patient can see the result at the same time in the double display.

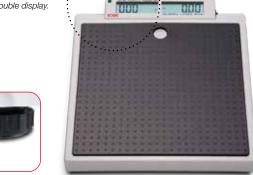



The flat design and slip-proof surface ensure easy access and a secure stance.



The large leveling bases on the bottom ensure stability on any type of flooring.



The two push buttons on the rear side of the scale can be operated with the foot.

seca 878

Technical Data

- Capacity: 200 kg
- Graduation: 100 g < 150 kg > 200 g
- Dimensions (WxHxD):
 321 x 61 x 362 mm
- Dimensions/platform (WxHxD):
 288 x 60 x 280 mm
- Weight: 4.1 kg
- Power supply: batteries
- Functions: mother/child-function, HOLD, automatic switch-off
- Optional: stadiometer seca 217
 combined with adapter element
 seca 437, carrying case seca 414
 (for combination with stadiometer and
 adapter element) and seca 415