seca mVSA 535



Spot-check Vital Signs Monitor

We make taking vital signs EASY



seca mVSA

Where simplicity meets accuracy

We know that the working environment for medical staff is complex. This is precisely why we developed the seca mVSA, a spot-check monitor that simplifies and quickens the routine task of capturing patient vitals (blood pressure, temperature, pulse rate, SpO_2) in less than 60 seconds.

- + EASY, quick and accurate
- + seca EasyQuick inflation blood pressure mode seca EQ BP[™]
- + Configurable seca EQ BP[™] average mode (BP30), sprint mode or single measurement
- + seca scale connectivity through seca 360° proximity
- + EMR-validated and certified
- + Incredible price-point

Component options

All measurements at a glance on the easy to operate 7" touchscreen

monitor.

+ World-class service program



Optional mobile stand seca 475.



Blood pressure cuffs (sizes XS-XXL)



Temperature probe (oral/axillary)



Temperature probe

(rectal)



Masimo RD SET® SpO₂ finger clip



Durable BP cuffs enhanced with metal plate and cuff stitching around the exterior.



Long lasting lithium-ion battery enables up to 150 measurements on one single charge.



Store up to 70,000 vital signs measurements in the large internal memory for later processing.

EASY, quick and accurate

The seca mVSA is easy-to-use and self-explanatory with little to no training required for nursing staff. Just turn it on, start the measurement and within 60 seconds blood pressure, temperature, pulse rate and SpO₂ are shown at a glance.





Record **blood pressure** with the unique inflation or deflation mode with a push of only one button in less than 20 seconds.



The **temperature** measurement starts automatically when the probe is used (oral/axillary or rectal).



The patient's **pulse rate** is always measured simultaneously while using blood pressure measurement or pulse oximeter.



Measure **blood oxygen levels** with a Masimo RD SET[®] sensor.

seca EQ BP™

True screening of hypertension with configurable seca EasyQuick inflation Blood Pressure mode, seca EQ BP[™]: take 6 readings automatically in 30 minutes with the seca EQ BP[™] average mode (BP30), or the average of 3 readings with the seca EQ BP[™] sprint mode. A single reading can be taken in under 20 seconds.



seca EQ BP[™] can be determined quickly and gently using the inflation measurement method.



Avoid the "white coat effect" and capture blood pressure accurately with configurable seca EQ BP[™] averaging or sprint mode.

Seamless scale connectivity

The seca mVSA can gather height and weight measurements from EMRvalidated seca scales and measuring stations, thereby improving and standardizing base parameter acquisition.





EMR-validated and certified


Eliminate transcription errors and increase patient safety with direct transmission of measurements into the electronic medical record (EMR) via HL7. Similar integration solutions are provided for a wide range of EMR systems:



Future safe with seca service

- + Advanced replacement service can immediately provide you with an equivalently configured unit to ensure minimal downtime
- + Priority access to latest seca software updates and upgrades
- + NEW: Optimize your investment with a convenient 2+5-year warranty extension
- + Outstanding Engineering Technicians dedicated to our nationwide customers. Operating times: 9-7 EST, 8-6 CST, 6-4 PST

Toll-free service hotline: +1 800 542 7322

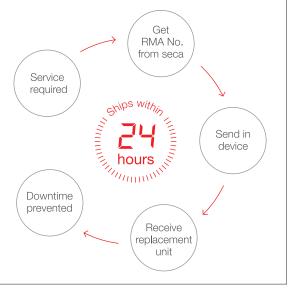

Configurations



Spot-check Vital Signs Monitor

- + NiBP 1 BP-Cuff size L
- + TEMP (oral/axillary)
- 1 Covidien Filac[™] 3000 Thermometer with oral/axillary probe 2 packages á 20 pcs
- 2 packages a 20 pcs Covidien Filac™ Probe Covers



Spot-check Vital Signs Monitor

+ NiBP 1 BP-Cuff size L

535 1020 131

+ TEMP (rectal) 1 Covidien Filac[™] 3000 Thermometer with rectal probe 2 packages á 20 pcs Covidien Filac[™] Probe Covers



Spot-check Vital Signs Monitor

+ NiBP 1 BP-Cuff size L



535 1200 131

Spot-check Vital Signs Monitor

+ NiBP 1 BP-Cuff size L

+ TEMP (oral/axillary)

1 Covidien Filac™ 3000 Thermometer with oral/axillary probe

2 packages á 20 pcs Covidien Filac™ Probe Covers

+ SpO₂

 Masimo RD SET[®] adult finger sensor
patient cable

535 1210 131



Spot-check Vital Signs Monitor

+ NiBP

1 BP-Cuff size L

+ TEMP (rectal)

- 1 Covidien Filac[™] 3000 Thermometer with rectal probe
- 2 packages á 20 pcs Covidien Filac™ Probe Covers
- + SpO₂
- Masimo RD SET[®] adult finger sensor
 patient cable

535 1220 131

535 1010 131

Accessories



Blood pressure cuffs

Size XS 5.1 – 7.9" / 13 – 20 cm	Size S 6.7 – 10.2" / 17 – 26 cm
490 0024 001	490 0025 001
Size M 9.4 – 12.6" / 24 – 32 cm	Size L 12.6 – 16.5" / 32 – 42 cm
490 0026 001	490 0027 001
Size XL 15.0 – 18.1" / 38 – 46 cm	Size XXL 16.5 – 21.3" / 42 – 54 cm
490 0028 001	490 0029 001



Extension tubes

490 0033 001

490 0034 001

Length: 51.2" / 130 cm

Length: 118.1" / 300 cm

Mobile stand Height adjustable: 27.4 – 46.6" / 69.5 – 118.5 cm

475 0535 009



Wall mounting kit Includes stand mount, cradle and GCX 5" slide-in plate

490 0017 001



Carrying bag

Dimensions (WxHxD): 16.1x12.2x11.6"/ 41x31x29.5 cm 432 0000 009



Filac[™] probe covers

Oral/axillary, rectal 2,000 pieces (100 x 20) (disposable)

490 0150 001

201353502022



seca north america | east 498 Wando Park Boulevard Mount Pleasant, SC 29464 USA

T +1 800 542 7322 F +1 888 705 7397 E info.us@seca.com

seca.com/EASY

